

Technical Data Sheet



1. Product


140.4 - Universal Indicator Strips pH 0 - 14, non-bleeding.

2. Instructions for use

- Immerse an indicator strip in the test solution for 1 second.
- Shake off excess solution.
- Wait 10-20 seconds.
- Compare the colour produced on the strip against the colour chart.

3. Measuring range

pH 0 - 1 - 2 - 3 - 4 - 5 - 6 - 7 - 8 - 9 - 10 - 11 - 12 - 13 - 14

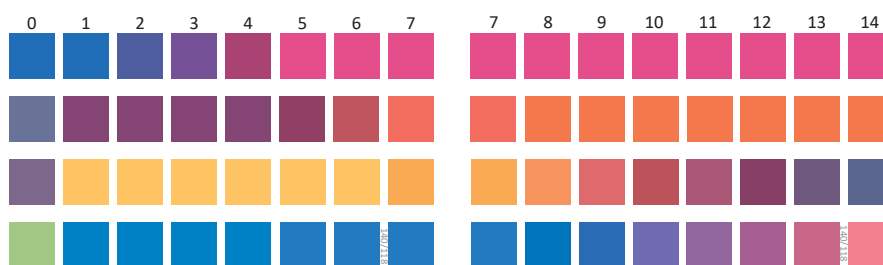

4. Application

pH indicator strip for testing of pH in solutions between pH 0 and pH 14.

5. Shelf Life

3 years from manufacturing date.

6. Contents

100 strips.

7. Weights & Dimensions

pH strip: 5m x 85mm

Plastic pack: 65mm(w) x 88mm(h) x 15mm(d)

8. Notes

Keep any unused strips in the original container. Do not store where they are liable to be exposed to acid or alkaline gases or vapours. Store in a cool dry place.

Johnson Test Papers Ltd
Hainge Park, Hainge Road
Oldbury, West Midlands
B69 2NU, United Kingdom

T: +44 (0)121 557 3883
F: +44 (0)121 557 8235
E: info@jtp.uk.com
W: www.johnsontestpapers.com

Registered in England
Company Reg No: 4112470
VAT No: GB 754 6913 04

